# Cellu M6° Integral [2] 1 Cellu M6° Endermolab [2] 1






















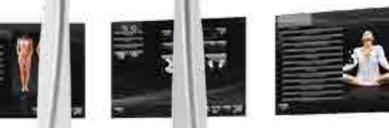













### LPG® GLOBAL PROTOCOL SYSTEM

### > Adjustable color touchscreen.

• More than 60 built-in protocols: Health / Sport / Physical Therapy.

# **ERGODRIVE HEAD**

- > New tissue massage +50%\* for more effective treatments.
- > Independent motorized rollers (LPG patents) 3 Roll rotations for precision treatments: Roll In, Roll Out, Roll Up.
- > Control screen and control buttons on the Ergodrive head and rotating base for easier use..

\*Compared to Keymodule technology

# **TR50 HEAD**

- > Specially designed for therapeutic use.
- > Independent motorized rollers (LPG patents).
- 3 adjustable tissue massage settings (deep to superficial) for targeted precision.
- **>** Control screen and control buttons on the TR50 head.

### **ERGOLIFT HEAD INNOVATION**

- > New LPG patent: integrated micro-motor [MPF: Motorized Pulsating Flaps] for two-dimensional Mécano-Stimulation™.
- > Tissue stimulation increased by 270%\*.
- **> Tissue constantly monitored** through the transparent chamber to control hyperemia.
- **> Easier access to narrow and sensitive areas** (cheek bones, underarms).
- \*Compared to LPG non-motorized FLAP technology

① Different depths of massage possible with sequential aspiration. ② Targeted cellular stimulation thanks to variable Flap pulsation speeds





# SCALABLE CONNECTIVITY

- > Plug & Care Technology for more versatility.
- **>** USB Port to **upload periodic LPG updates** (new protocols).



Advanced aesthetic and therapeutic technologies





Weight: 168 lbs
Dimensions (Lxwxh): 31x27x66 inches
Frequency: 2 - 16 Hz
Pressure: 690 mbar
Output: 26.15yd³/h
Sound level: 56 dBA
Electrical characteristics:
100-240 V / 50-60 Hz / 625-650 W
Color touchscreen: 10"4
Energy saving mode: yes
Audio: yes
USB Port: yes

Endermologie® is an exclusive technology of cellular stimulation using Mécano-Stimulation™.

Created in 1986 by LPG®, it non-aggressively treats connective tissue for beauty or therapeutic purposes.

Its effectiveness has been backed by more than

110 scientific studies.

Cellu M6° Endermolab [2] ii

Cellu M6° Integral [2] 🗓



# **ERGOLIFT HEAD**

**FLAP 10:** narrow areas around eyes and lips, hands and fingers.

**FLAP 20:** large areas of thin and sensitive tissue.



# **AUXILIARY HEADS**

To reach the most narrow and sensitive areas (lesions, scars, etc.).



### **MICRO-HEADS AND MICRO-NOZZLES**

For micro-targeted use (fingers for example).




# **ERGODRIVE HEAD**

Driving seat of the system equipped with KM80 for deep action on large areas.



### **TR50**

For therapeutic use: fibrosis, edema, inflammation, pre- and post-liposuction, etc.

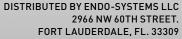

Toll Free: [866] 374-9401 Local: [954] 636-3000 Fax: [954 984-8836 www.endosystems.us

